

eLINE[®] 0.1-5 µl


with Super Pipetting feature

Now you can dispense the smallest volume easily and precisely! eLINE 0.1 µl - 5 µl electronic pipette dispenses the very smallest volumes from the air - no need to touch the receiving vessel wall or dispense under the liquid surface


- Contamination free pipetting
- Improved accuracy in small volumes
- Less tip consumption – savings in tip costs
- Ideal in PCR, microplate pipetting

Extremely rapid piston movement allows complete emptying of the tip!

Precise dispensing in small volumes

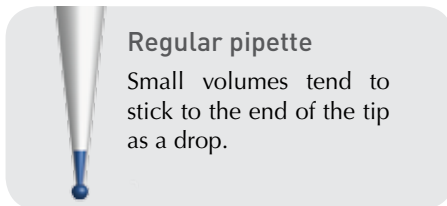

Regular pipette
Small volumes tend to stick to the end of the tip as a drop.

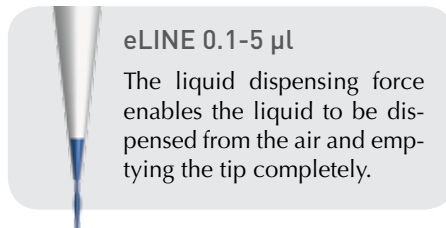

eLINE 0.1-5 µl
The liquid dispensing force enables the liquid to be dispensed from the air and emptying the tip completely.



WATCH VIDEO ON
www.biohit.com

eLINE 0.1-5 µl operation modes:

- Super Pipetting (SP)
- Pipetting (P)
- Reverse Pipetting (rP)
- Pipetting with Cycle Counter (Pcount)
- Multiple Dispensing (d)
- Diluting (dd)
- Diluting with Mixing (dd+Mixing)
- Automatic Multiple Dispensing (Ad)
- Multi-Aspirating (SA)
- 6 memory locations

Performance specifications and ordering details

Cat. No.	Channels	Volume Range [µl]	Increment [µl]	Test Volume	Inaccuracy %	Imprecision %	Compatible Biohit Tip
73001X	1	0.1-5 µl	0.05	0.5	9.50	8.00	10 µl
				2.5	1.70	1.20	
				5.0	1.20	0.70	

0=without AC-adaptor1, 1=with universal AC-adaptor

eLINE Charging Stand

Cat. No.	Item
730981	eLINE Charging Stand*
730991	eLINE Charging Carousel*

* Incl. universal AC-adaptor

